







JOYCEeasy

...ingenious simplicity.

Nasal or Full Face. Simply easy.

It could not be any simpler. Everything about our ventilation mask JOYCEeasy is designed for easy use – from simple mask adjustments directly on the patient's face to intuitive handling. Because the mask fits so well, JOYCEeasy also encourages better therapy compliance by minimizing leaks and pressure points. The quiet, diffuse exhalation system ensures a peaceful bedroom atmosphere.

Ingenious simplicity everywhere:

- Unlimited freedom of movement provided by the ball-and-socket joint
- Quiet and diffuse exhalation system 2
- Mask can be put on and taken off with just one hand
- Color-coded hook-and-loop tabs 3
- Extremely robust plastic
- · Quick-release cord (only Full Face)

JOYCEeasy next FF JOYCEeasy Full Face

... the perfect duo.

The right fit for everyone.

We heard you.


And we reworked JOYCEeasy Full Face to satisfy your practical needs.

JOYCEeasy next FF now features a completely redesigned large mask cushion that has a comfortable and firm fit even at high pressures.

The headgear too remains securely in place at high pressures. Thanks to the right material mix and ergonomic design, the headgear clips can be easily loosened from the mask.

Your benefit – The right mask for everyone:

JOYCEeasy next FF = normal to obese JOYCEeasy Full Face = slender to gaunt

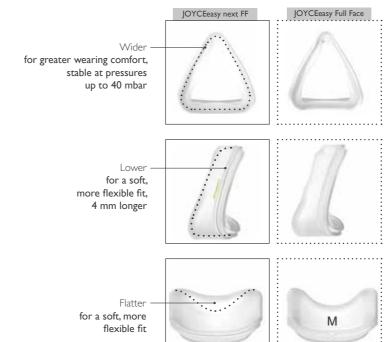

JOYCEeasy nasal

Thanks to your feedback, we have made more improvements to the JOYCEeasy mask cushion.

The wider and softer cushion sits lower across the bridge of the nose and easily accommodates larger noses comfortably. The flatter and more flexible upper lip area lies more gently on the face.

Now even more patients can benefit from the enlarged mask cushion.

Don't worry! The size of the cushion has not been changed in Size S, which still fits small faces.

Patient Interface / JOYCEeasy

Patient Seasy



Ball-and-socket set, vented Full Face (JOYCEeasy next FF) anti-asphyxia valve, rotating sleeve WM 15944

Ball-and-socket set, vented Full Face (JOYCEeasy Full Face) retaining ring, anti-asphyxia valve, rotating sleeve WM 15897



Endoscope adapter set JOYCEeasy next FF NV/ JOYCEeasy Full Face NV retaining ring NV, endoscope elbow NV WM 15864

Elbow set JOYCEeasy next FF NV / JOYCEeasy Full Face NV retaining ring NV, elbow NV WM 25269

Replacement parts



WEINMANNstrap for JOYCEeasy Full Face WM 25606

Mask cushion

Size S WM 26803 (Nasal) WM 26853 (Full Face)





_	
	Forehead cushion silicone WM 26801 (Nasal and Full Face)

Product name	Size S	Size M	Size L
JOYCEeasy Nasal	WM 25510	WM 25520	WM 25530
JOYCEeasy Full Face	WM 25910	WM 25920	WM 25930
JOYCEeasy next FF	WM 25109	WM 25119	WM 25129
JOYCEeasy next FF NV	WM 25159	WM 25169	WM 25179

For more information on our therapy solutions, accessories and mask systems please visit loewensteinmedical.de

C€ 0197





I E-Centre, Easthampstead Road

Löwenstein Medical UK

RG12 | NF Bracknell

United Kingdom

T: +44 1344 830023

info@loewensteinmedical.uk

www.loewensteinmedical.uk

Löwenstein Medical Technology Kronsaalsweg 40, 22525 Hamburg Germany T: +49 40 54702-0 F: +49 40 54702-461 info@loewensteinmedical.de www.loewensteinmedical.de



JOYCEeasy

Easy to Fit, Easy to Use.



