


## MOVING ELEGANCE

NeuroSlider® Microcatheter DLC

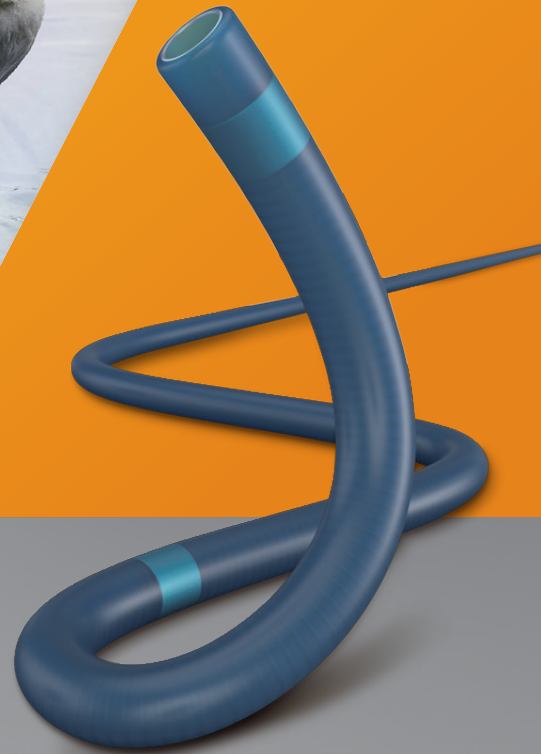

- Superior torqueability and pushability
- Smooth and safe device delivery
- Longlasting tip shape retention

## ADVANCE.


- Dual layer hydrophilic coating ensures outstanding lubricity and durability.
- Braiding / coiling reinforcement induces superior torqueability and significant reduction of ovalisation and elongation.

## NAVIGATE.

- Shapeable tip with lasting shape retention allows excellent distal navigation even in tortuous anatomies.
- Multi polymer construction consisting of different flexibility zones with smooth transitions from maximum stability at the hub to maximum flexibility at the tip permits precise and effective navigation.

## DELIVER.

- Inner PTFE liner minimises friction and allows a controlled and safe delivery of therapeutic and diagnostic agents.
- Advanced hub design with a transparent window results in a precise device transfer into the hub.

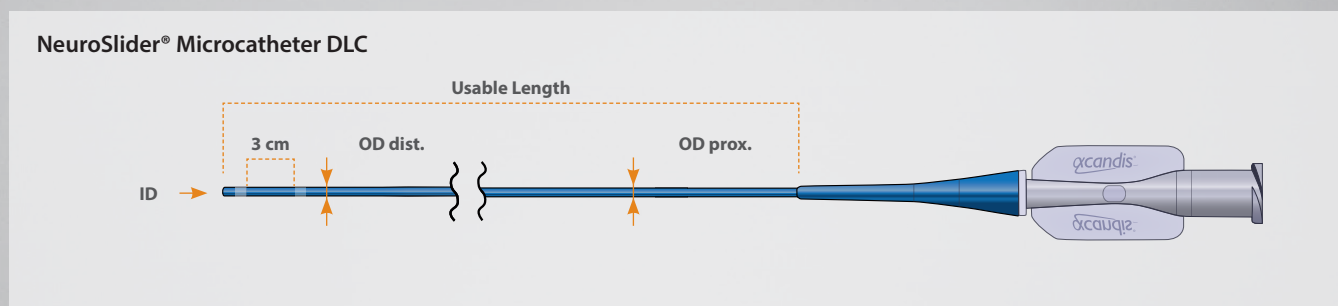

## ORDERING INFORMATION

Product Name	Reference Number	ID (Inch)	OD dist. / prox. (French)	Usable Length (cm)	Tip Shape	Tip Marker
NeuroSlider® 17	01-000272	0.0165	1.9 / 2.1	155	Straight (shapeable)	2
NeuroSlider® 17 DLC	01-000282	0.0165	1.9 / 2.3	155	Straight (shapeable)	2
	01-000283	0.0165	1.9 / 2.3	160	Straight (shapeable)	2
	01-000284	0.0165	1.9 / 2.3	167	Straight (shapeable)	2
NeuroSlider® 21	01-000273	0.021	2.4 / 2.5	155	Straight (shapeable)	2
NeuroSlider® 21 DLC	01-000292	0.021	2.2 / 2.6	155	Straight (shapeable)	2
	01-000293	0.021	2.2 / 2.6	160	Straight (shapeable)	2
	01-000294	0.021	2.2 / 2.6	167	Straight (shapeable)	2
NeuroSlider® 27 (DLC)	01-000274	0.027	3.0 / 3.1	155	Straight (shapeable)	1

All changes or modifications, may they be technical or other, or changes in the availability of products are expressly reserved.

Distributed by:

**acandis®**  
ENGINEERING STROKE SOLUTIONS

CE 0297

**ACANDIS GmbH**  
Theodor-Fahrner-Str. 6  
75177 Pforzheim  
Germany

Tel: +49 7231 155 00 0  
Fax: +49 7231 155 00 129  
E-Mail: [info@acandis.com](mailto:info@acandis.com)  
[www.acandis.com](http://www.acandis.com)