







HealFix® Ti Suture Anchors Family

Various options in anchor design, diameter, suture size/color, with/without needles



reddot winner 2020

The Red Dot Award:Product Design only awards this sought-after seal of quality to products that feature an outstanding design. In 2020, HealFix Suture Anchor won the Red Dot award in the category "Medical Devices and Technology".

-  1.3mm with needles
2.0mm with needles
2.4mm with needles
-  2.0mm with needles
2.0mm
-  2.8mm with needles
2.8mm
-  3.5mm with needles
3.5mm
-  5.0mm with needles
5.0mm
6.5mm with needles
6.5mm
-  4.5mm Plus
5.5mm Plus
6.5mm Plus

rejoin 

Hangzhou Rejoin Mastin Medical Device Co., Ltd

Headquarter Office: 22 Xinyan Rd, Yuhang Economic Development Zone, Hangzhou, Zhejiang, 311100, P.R. of China
Beijing Branch: Room 1605, Jinghui building, No. 118 Jianguo Road, Chaoyang District, Beijing, P.R. of China
Shanghai Branch: 309A, L1-B, Hongqiao World Center, 1588 Zhuguang Rd, Qingpu District, Shanghai, P.R. of China

contact@rejoin.com @rejoinmedical www.rejoin-medical.com

Internal information is for professional reference only

rejoin 




HealFix® Anchor
Titanium Fixation System

Wide Range of Applications

More options for sizes, threads & sutures

Pre-loaded with WilSuture® Ultra sutures

Ergonomic handle design provides surgeons comfort and ease of use

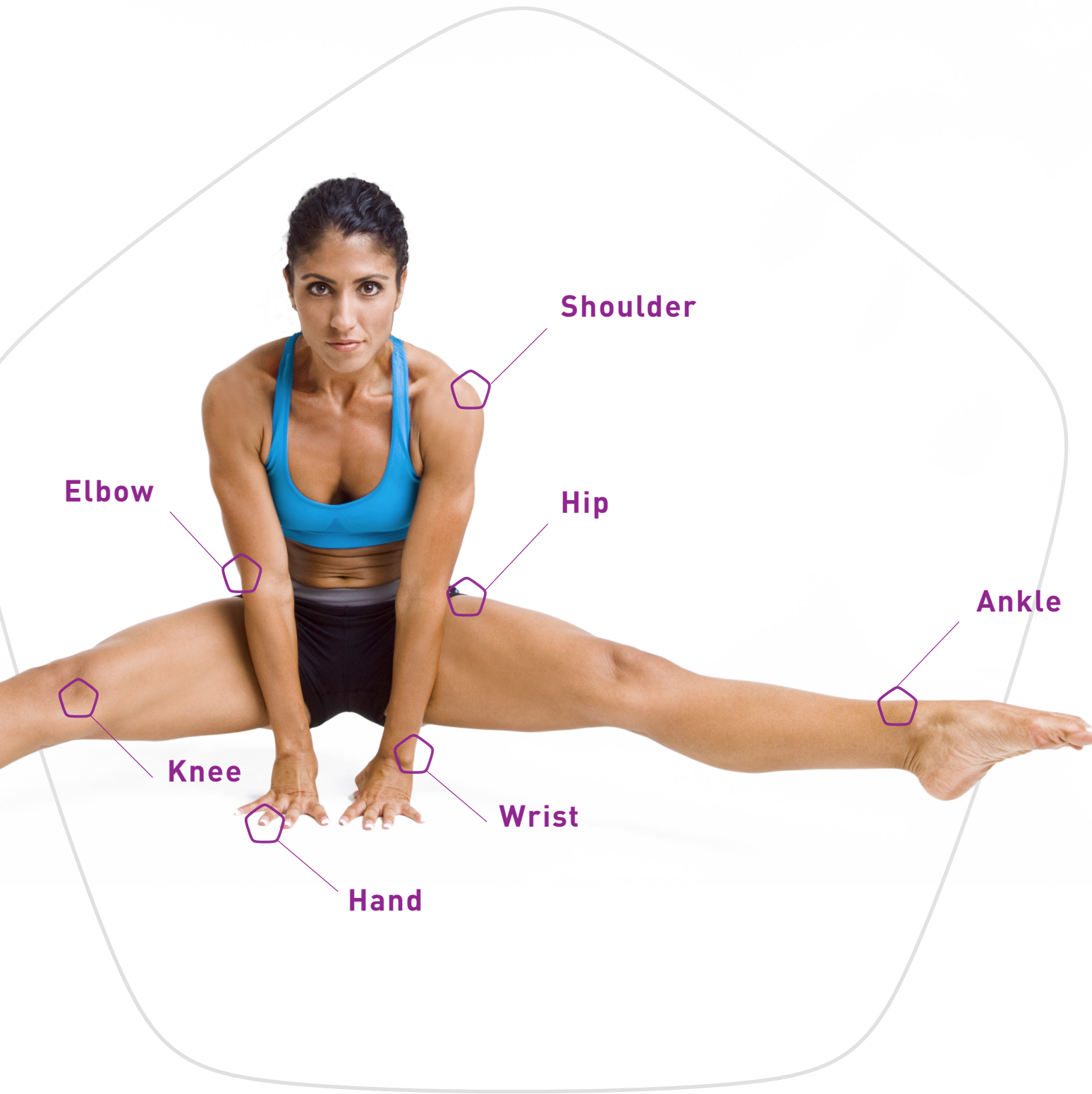

HealFix® Suture Anchor, Ti

- High-Low double thread design enhances the pull-out strength
- Two independent eyelets provide optimal fixation of both sutures & eliminate suture twining
- Sharpened tip facilitates easy insertion



HealFix®

HealFix® Plus Suture Anchor, Ti

- Fully Helix transitional thread design increases contact areas with cortical bone for high pull-out strength
- Built-in suture hang-up hole eliminates Cheese-Wire effect between sutures and bone after screw implantation
- Self-punching design makes implantation easier




HealFix® Plus

HealFix® Mini Suture Anchor, Ti

- NiTi arcs design provides subcortical purchase with outstanding pull-out strength
- Smaller size allows comprehensive fixation and restoration of natural footprint
- 1.3mm minimum size provides optimized soft tissue fixation for DEX applications



HealFix® Mini

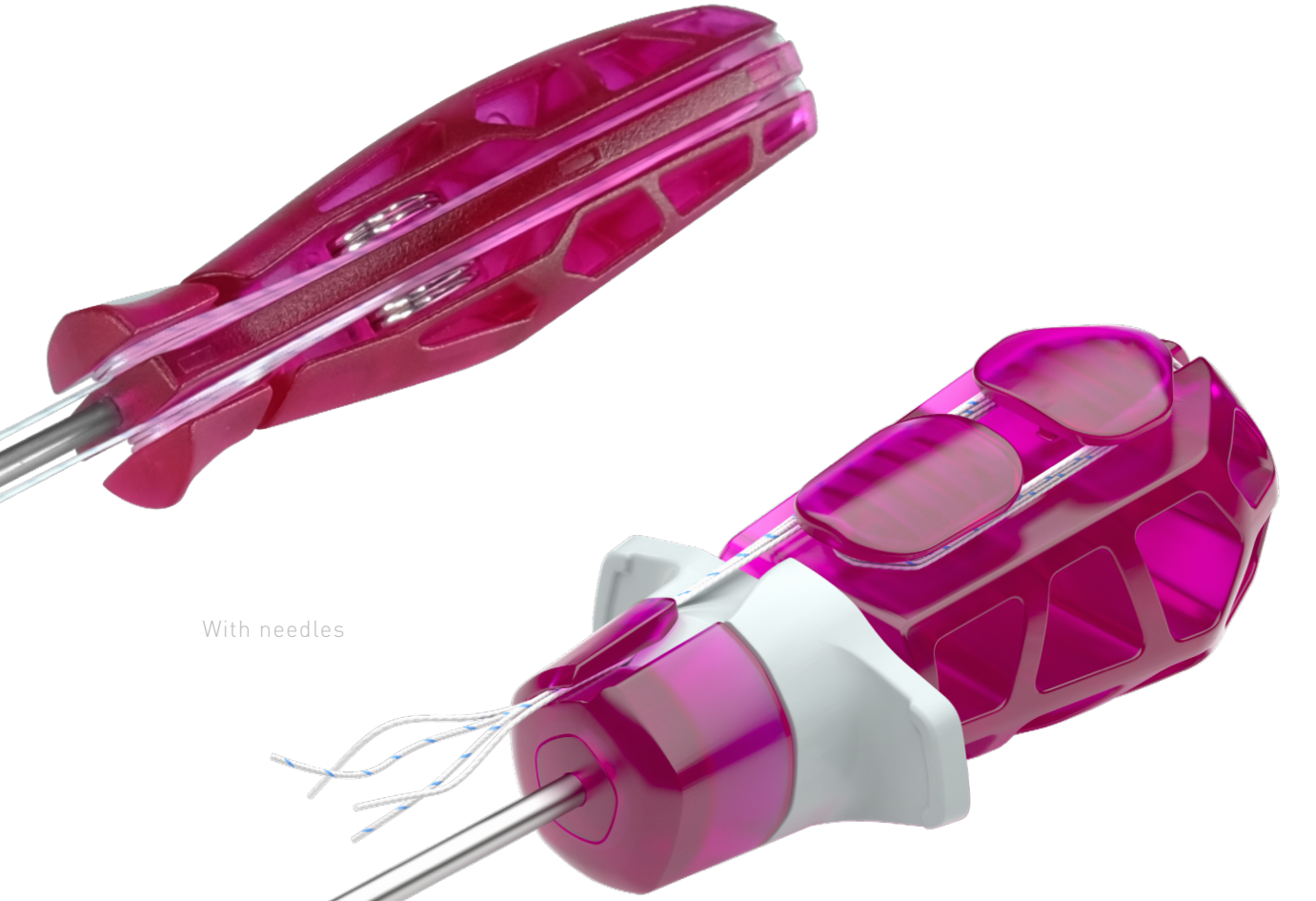

With needles

Without needles