



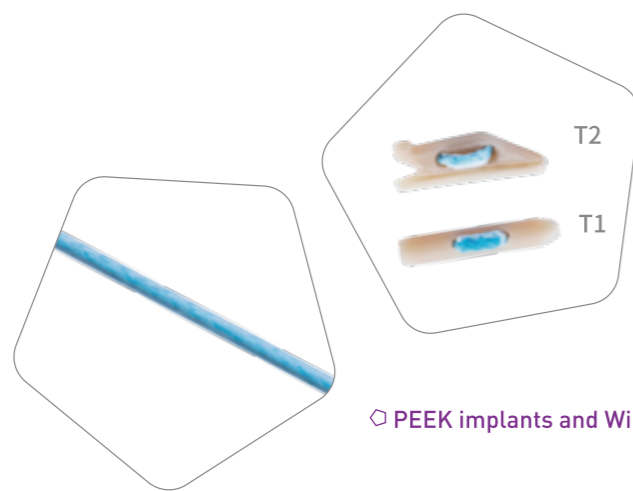


RapidLink® Meniscal Repair System


RapidLink® Meniscal Repair System consist of two PEEK implants and #2-0 UHMWPE Wilsuture, which provide a strong & simple minimally invasive solution for meniscal repair

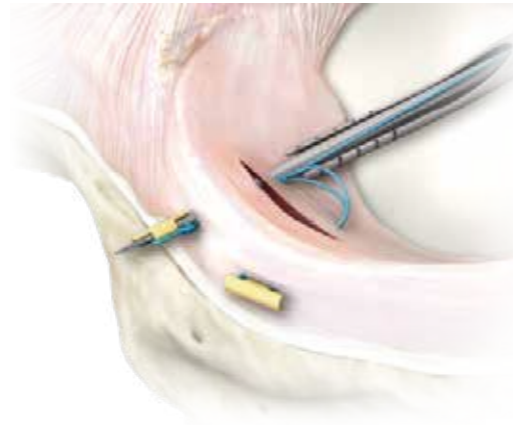


◊ PEEK implants and Wilsuture Ultra #2-0

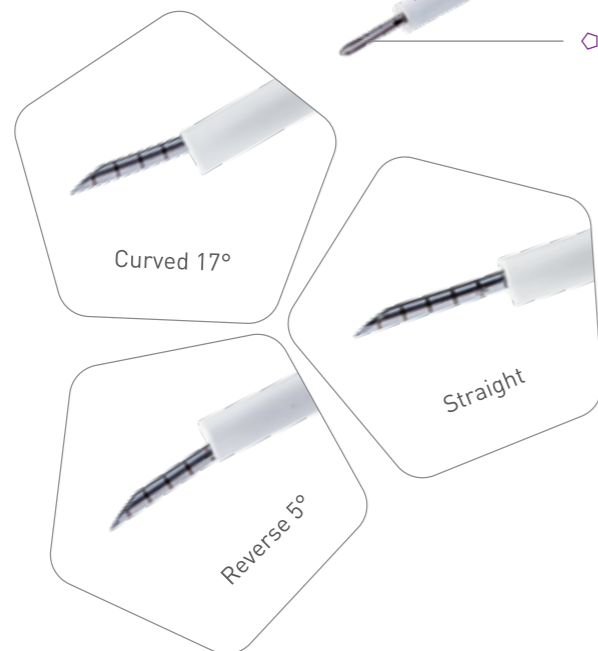

◊ Various angles, depth laser mark

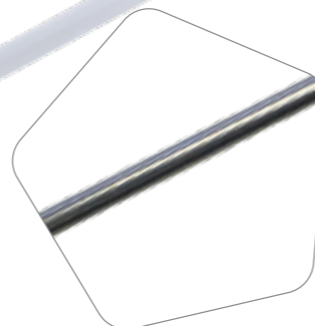

◊ Reinforced needle rod

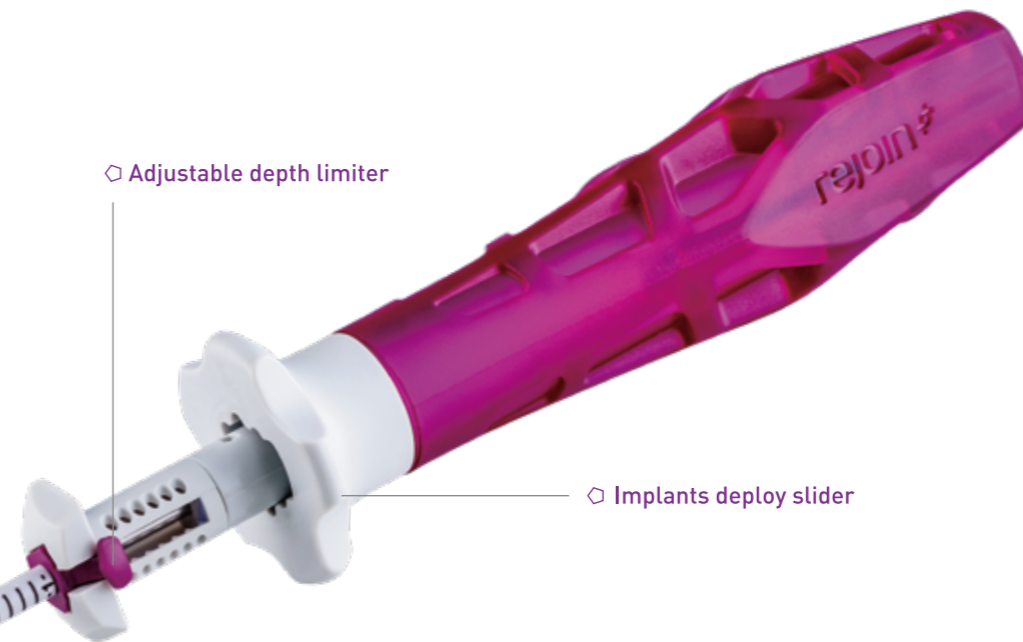


◊ Adjustable depth limiter

◊ Implants deploy slider

RapidLink® Plus All-inside Repair System

- ◊ **Minimally Invasive**
Smaller implants and Wilsuture Ultra #2-0 minimize disruption to the meniscus
- ◊ **Convenience**
360° active deploy implants at any hand position.
Auditory confirmation assures that implants are deployed
- ◊ **Detail**
Adjustable depth limiter eliminates the risk of early deployment and damaging nearby neurovascular structures
- ◊ **High Efficiency**
Reinforced needle rod penetrates with ease

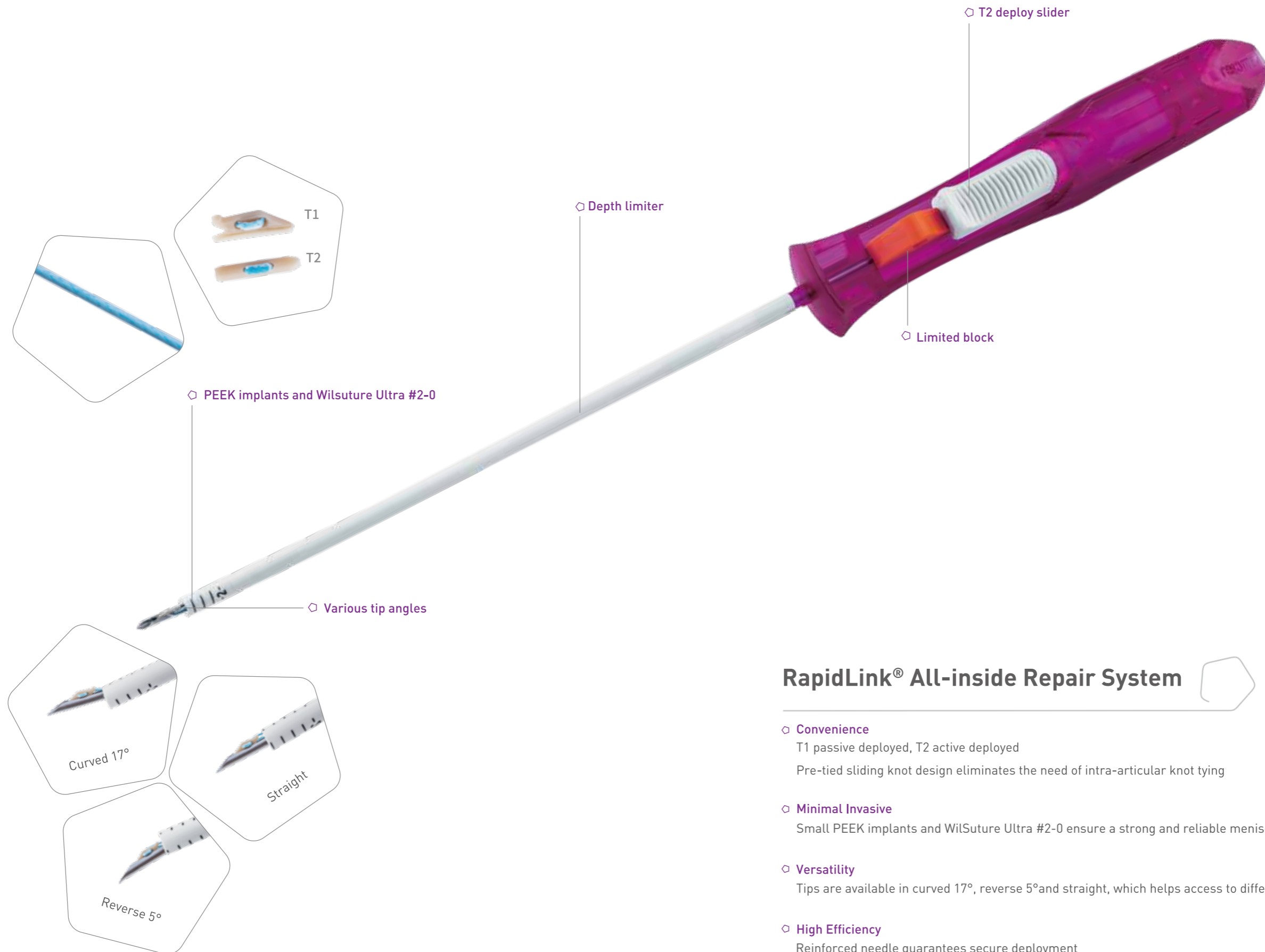

RapidLink® All-inside Repair System

- ◇ **Convenience**
 - T1 passive deployed, T2 active deployed
 - Pre-tied sliding knot design eliminates the need of intra-articular knot tying
- ◇ **Minimal Invasive**
 - Small PEEK implants and Wilsuture Ultra #2-0 ensure a strong and reliable meniscal repair
- ◇ **Versatility**
 - Tips are available in curved 17°, reverse 5° and straight, which helps access to different zone of meniscus.
- ◇ **High Efficiency**
 - Reinforced needle guarantees secure deployment



TwinAccord® Knot Pusher & Suture Cutter

- ◊ Arc-shaped edge designed to cut suture efficiently
- ◊ Titanium Nitride coating ensures a sharp and durable cutting edge
- ◊ Suitable for #0 & #2-0 sutures



TwinAccord® Meniscal Tool

- ◊ Unique design suitable for meniscus tear edges refreshing
- ◊ Multi-helical tooth design to freshen the edge of the meniscal tear
- ◊ Ergonomically designed to access all areas of the meniscus



Code	Description	Code	Description
800100100	RapidLink Plus Meniscal All-inside Repair System, Straight	800101200	Meniscal Depth Probe
800100200	RapidLink Plus Meniscal All-inside Repair System, Curved	800101300	45° Diamond Rasp
800100300	RapidLink Plus Meniscal All-inside Repair System, Reverse Curved	800101400	90° Diamond Rasp
800100400	RapidLink Meniscal All-inside Repair System, Straight	800101100	Slotted Cannula
800100500	RapidLink Meniscal All-inside Repair System, Curved		
800100600	RapidLink Meniscal All-inside Repair System, Reverse Curved		
810612100	Slotted Cannula, Lite		
810600200	Knot Pusher/Suture Cutter, 2-0# & 0#		
810610200	Knot Pusher/Suture Cutter, 2-0# & 0#		